



Quick Start: Mounting a ProEM+:1600 Series (3.60 3-hole Slotted Flange) to an Acton Series Spectrograph

Tools:

5/32" Hex wrench

3/32" Hex wrench

5/16" Open end wrench (provided in Spectrometer Kit bag)

(3) 10-32x 1/2" hex head stainless steel screws (provided in Spectrometer Kit bag)



Procedure:

1. Remove shipping plate.



2. Loosen setscrew with 3/32" hex wrench.



3. Loosen setscrew.



4. Gently rotate and pull.







5. Finger tighten hex head screws.
Leave about ¼" thread exposed.



6. Mount tube to camera, open shutter,
and align baffle with CCD array.

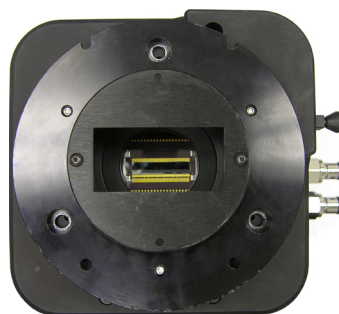

7. Tighten mounting screws.

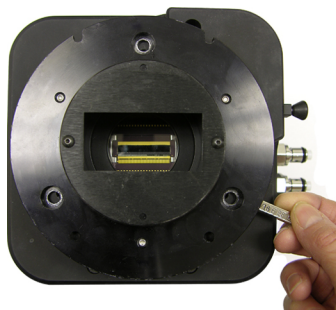

8. Gently rotate while inserting.



9. Tighten setscrew.



10. Tighten setscrew.

