



15 Discovery Way, Acton, MA 01720

Phone: (978)263-3584, **Fax:** (978)263-5086

Web Site: www.piacton.com

Supplemental Instructions
CCD Camera Adapter Plate to an Acton Series
Monochromator/Spectrograph
for C-Mount CCD Camera

Revision History

<i>Revision</i>	<i>Reason For Change</i>	<i>ECN</i>	<i>Date Issued</i>
A	Add "Remove outer baffle from sliding tube when using with SP-215X or SP-275X".	4157	11-10-06

Per customer specification, this Acton Series monochromator/spectrograph is supplied with an Adapter Plate to mount the above specified CCD with a focal plane distance of 17.5 mm to its mounting flange. To mount the CCD to the instrument's exit port the following procedure is recommended:


1. Remove the slide tube with the attached Adapter Plate from the monochromator/spectrograph following the instructions in the supplied manual.
2. Become familiar with the sliding tube and adapter assembly noting the labels "TOP". Remove the Adapter Plate from the sliding tube by removing the three (3) #30-32 button head screws (1/8" Allen wrench).
3. Remove the tape from the Adapter Plate/Clamp Screw.
4. Attach the CCD to the Adapter Plate with the clamp screw. Assure the orientation is as shown making sure the label "TOP" is aligned with the top of the CCD and securely tighten the clamp screw.
5. Reattach the adapter plate/camera assembly to the sliding tube with the labels "TOP" at the top.
6. Follow the instructions in the supplied manual for installing the sliding tube into the end exit port of the instrument. Note: The gap between the front of the monochromator/spectrograph and the sliding tube flange should be set at approximately 1.3 mm for best focus.
7. Remove outer baffle from sliding tube when using with an SP-215X or SP-275X.

CAUTION!

The CCD must be supported so as not to apply any undue stress or strain to the instrument.

THIS DRAWING REMAINS THE PROPERTY OF
PRINCETON INSTRUMENTS AND
CONTAINS PROPRIETARY INFORMATION WHICH
SHALL NOT BE REPRODUCED OR DISCLOSED
WITHOUT WRITTEN AUTHORIZATION.

REV	DESCRIPTION	DATE
1	.153 WAS .201	05-16-01

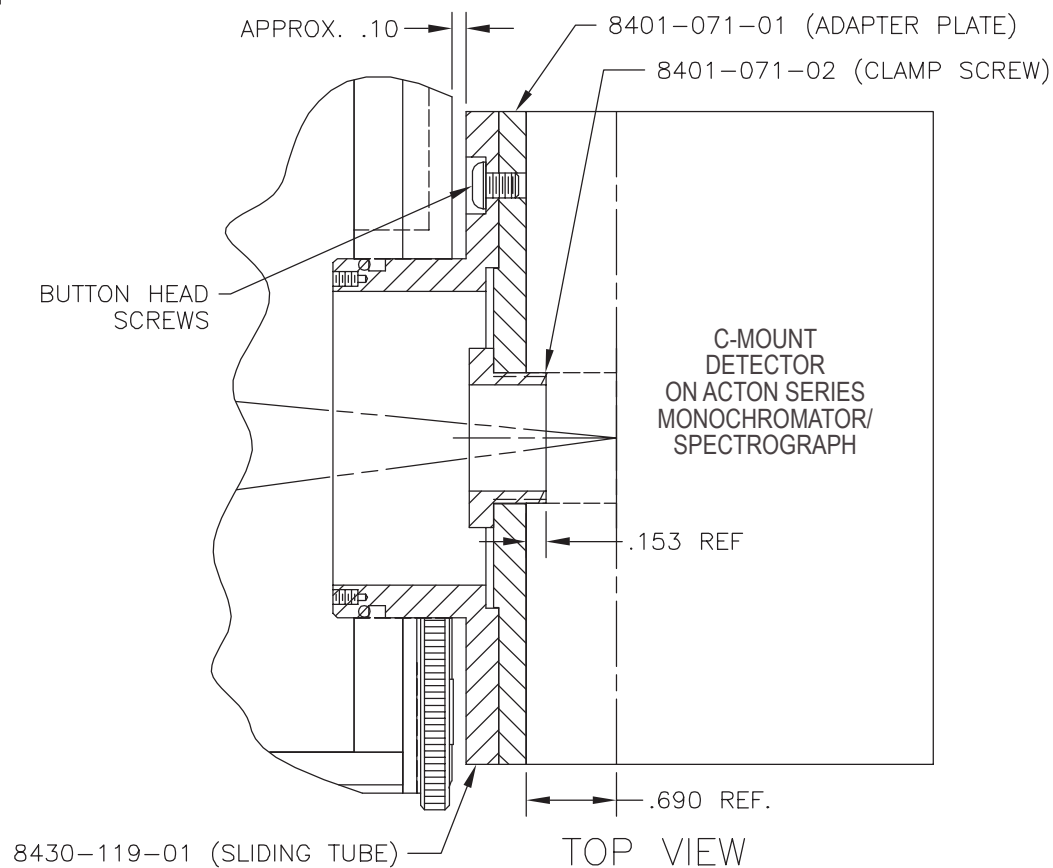

UNLESS OTHERWISE SPECIFIED		MATERIAL		FINISH	
DIMENSIONS ARE IN INCHES TOLERANCE ON FRACTIONS DECIMALS ANGLES $\pm 1/64$.XX $\pm .010$ ± 15 MIN .XXX $\pm .005$ MACHINE SURFACE 32✓ MAX REMOVE ALL BURRS & BREAK SHARP EDGES R.010 MAX		SCALE: FULL		JO#34539	
		DATE: 10-2-98		DRAWN BY: HCL	
		NEXT ASS'Y:		QTY: 1	
				15 DISCOVERY WAY ACTON, MA 01720	
DWG. TITLE		SIZE	DWG NO.	REV.	
C-MOUNT		B	8401-071-03	1	