



15 Discovery Way, Acton, MA 01720

Phone: (978)263-3584, **Fax:** (978)263-5086



Web Site: www.piacton.com

**Operating Instructions
for
SC-447
Sample Chamber**



**SC-447 Sample Chamber
mounted on the exit slit of the
SP-150 scanning monochromator.**

The SC-447 is compatible with all Acton Series monochromators and spectrographs (mounts to the standard slit assembly).

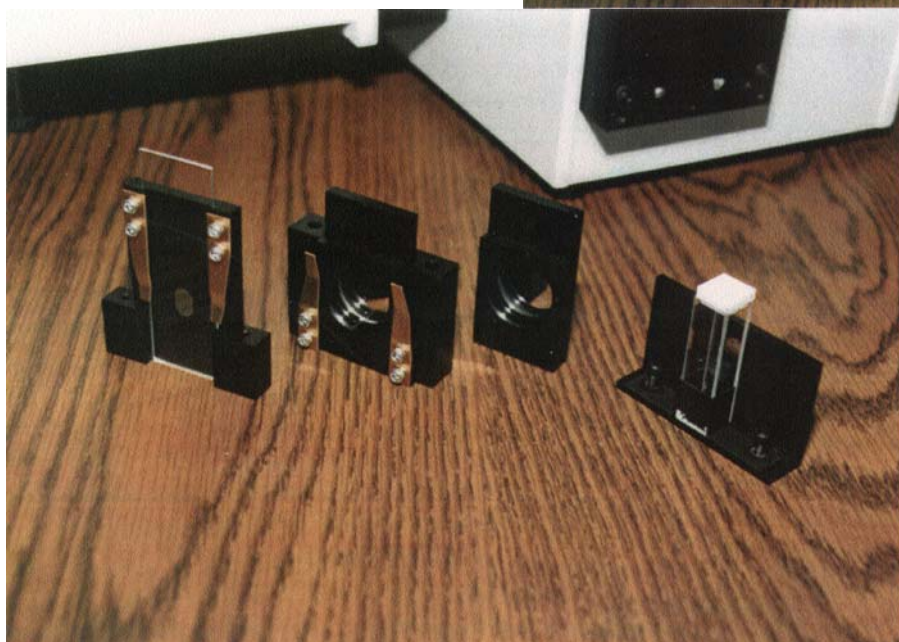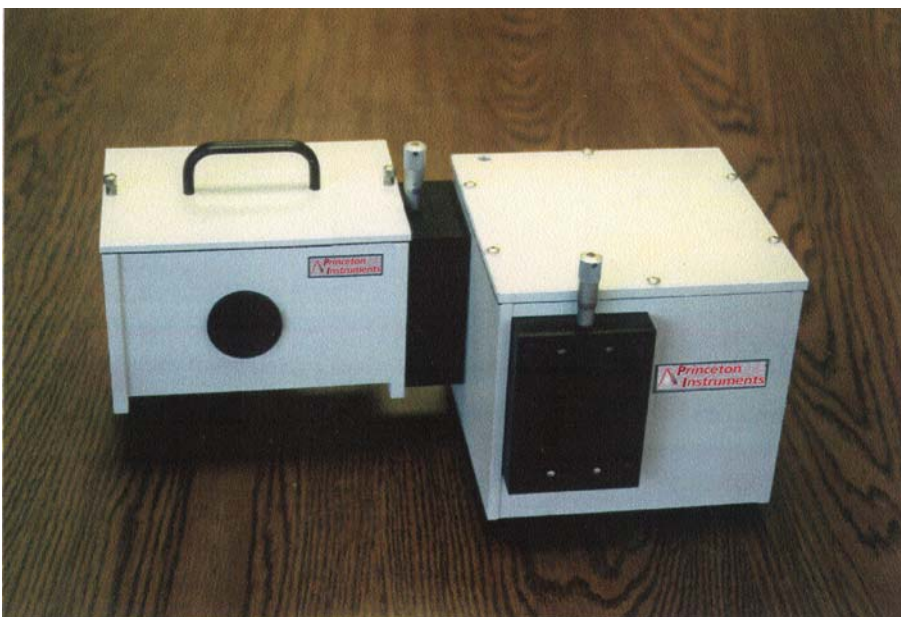

**Standard Sample Holders
(shown from left to right)**

447-3 Fixed Position Microscope
Slide Holder

447-2 Fixed Position Mount with
Removable Sample Holder for
1.0" (25.4 mm) Diameter Samples


447-2A Additional 1.0" (25.4 mm)
Sample Holder


447-1 General Purpose Holder,
shown with standard 110 mm
path length cuvette.



**SC-447 Sample Chamber with
the 447-1 General Purpose
Holder set up for standard
10 mm path length cuvettes.**

Notice the optical configuration (from left to right) showing the monochromator exit slit, order sorting filter holder, first quartz lens, 447-1 general purpose holder mounted on a rotary carousel, a second quartz lens, and the sample chamber output port with tapped holes for mounting detector assemblies.

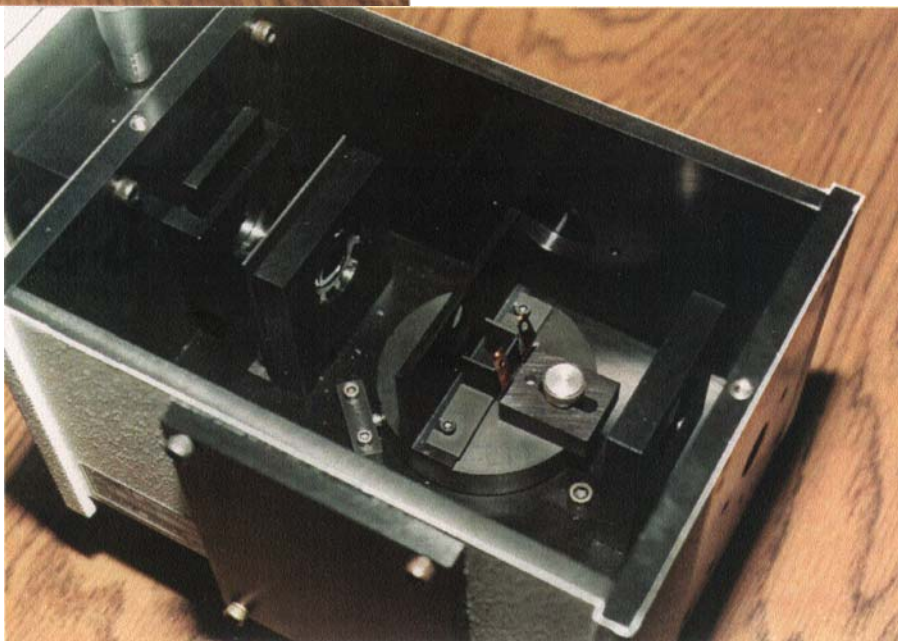


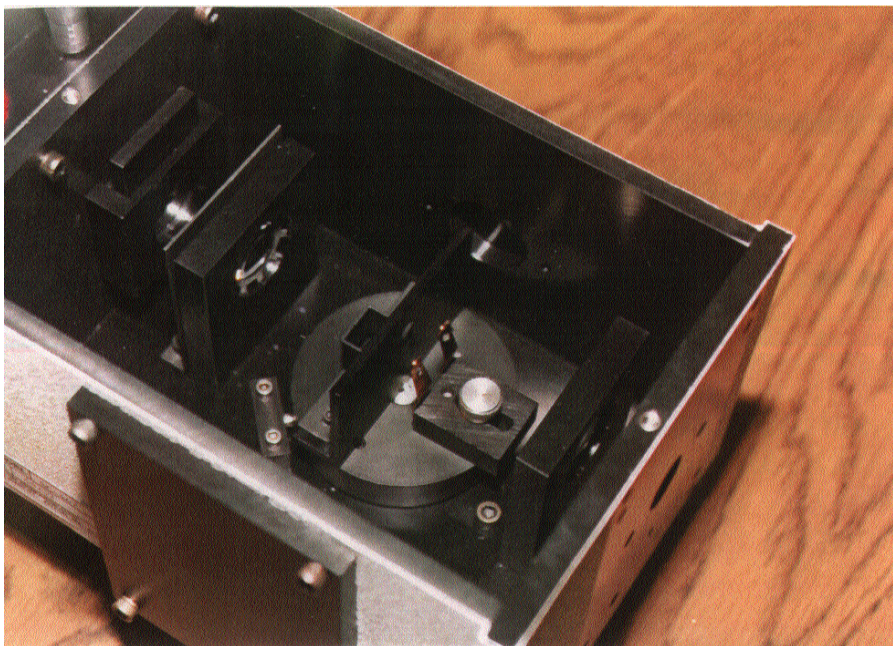


Inside view of the SC-447 showing the 447-1 general purpose sample holder turned around to accept solid samples up to 50x50x10mm.

Notice the "spring clip" style holders mounted on the adjustable sliding block to hold solid samples securely in position.



SC-447 with the 447-2 fixed position mount and removable sample holder for 1.0" (25.4mm) diameter samples.

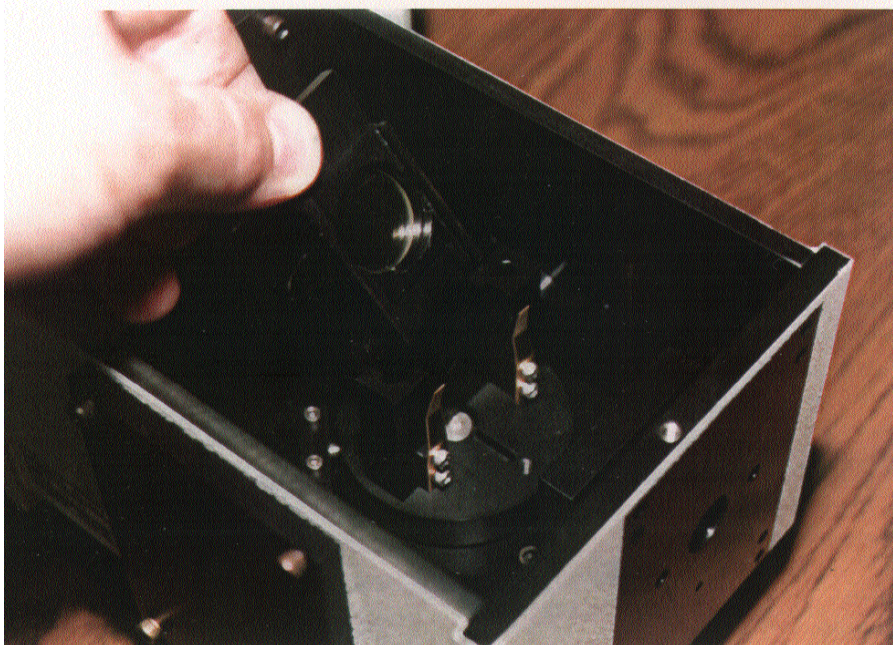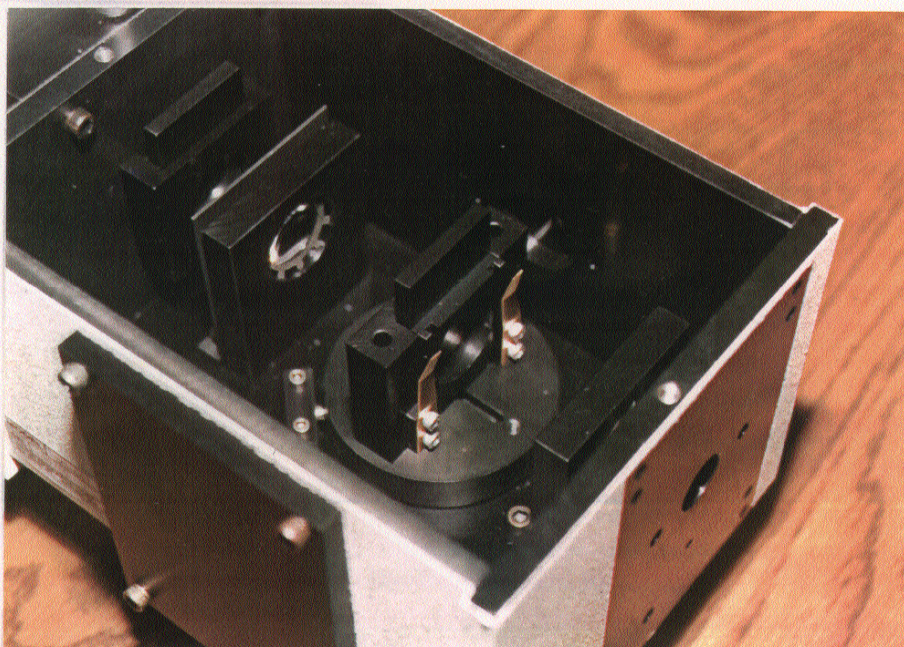


The 447-2A sample holder being removed from the fixed position mount.

The sample holder is designed for easy insertion and removal, for fast, convenient sample changing.


MODEL SC-447 SAMPLE CHAMBER

DETECTOR / SOURCE POSITION AND
ACCESSORY HOLE LOCATIONS
FOR PORTS 2 & 3

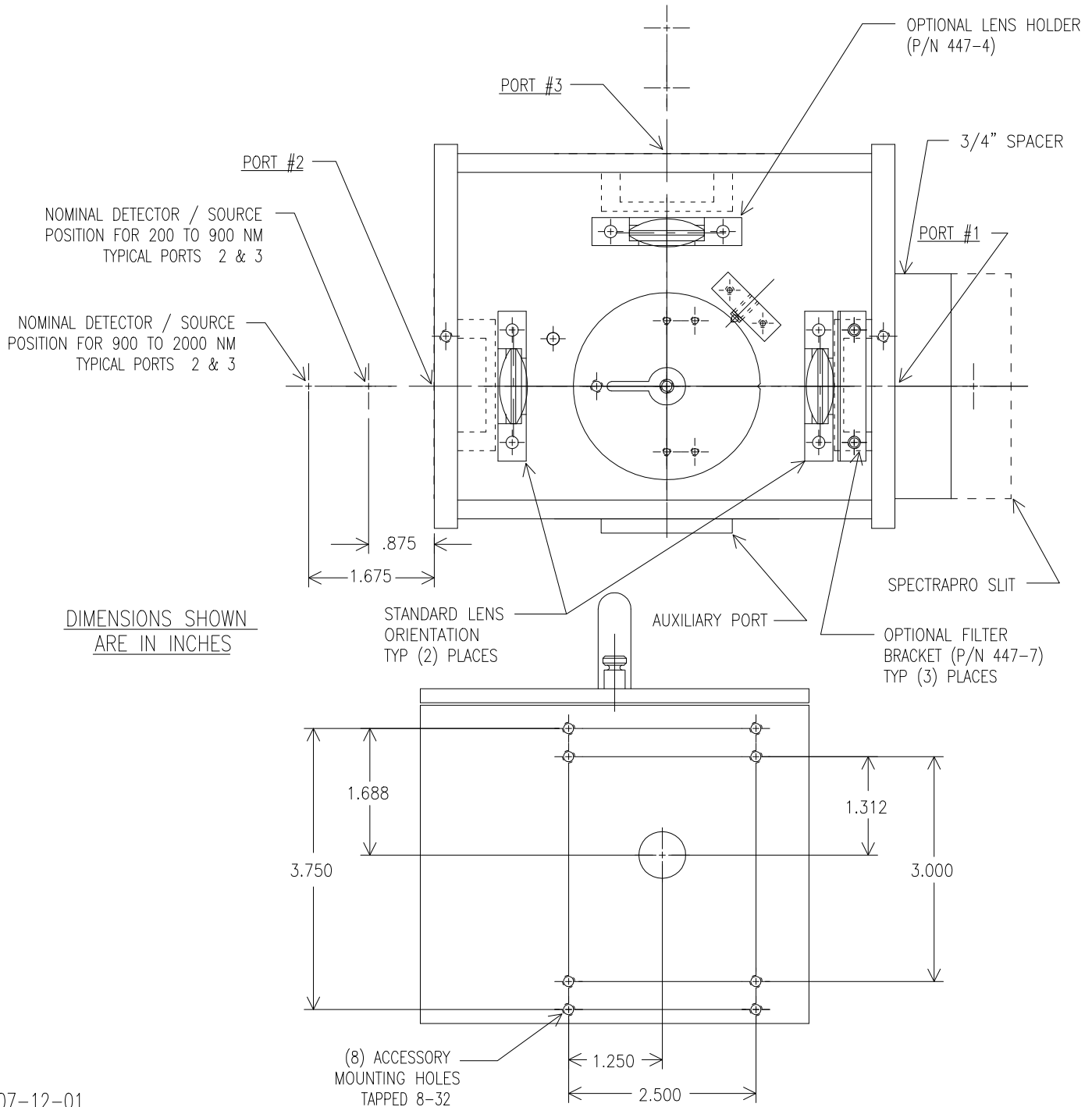

MODEL 447-1 GENERAL PURPOSE SAMPLE HOLDER
 CUVETTE HOLDER AND APERTURE BRACKET POSITIONS
 IN MODEL SC-447 SAMPLE CHAMBER

SOLID SAMPLE
 CONFIGURATION

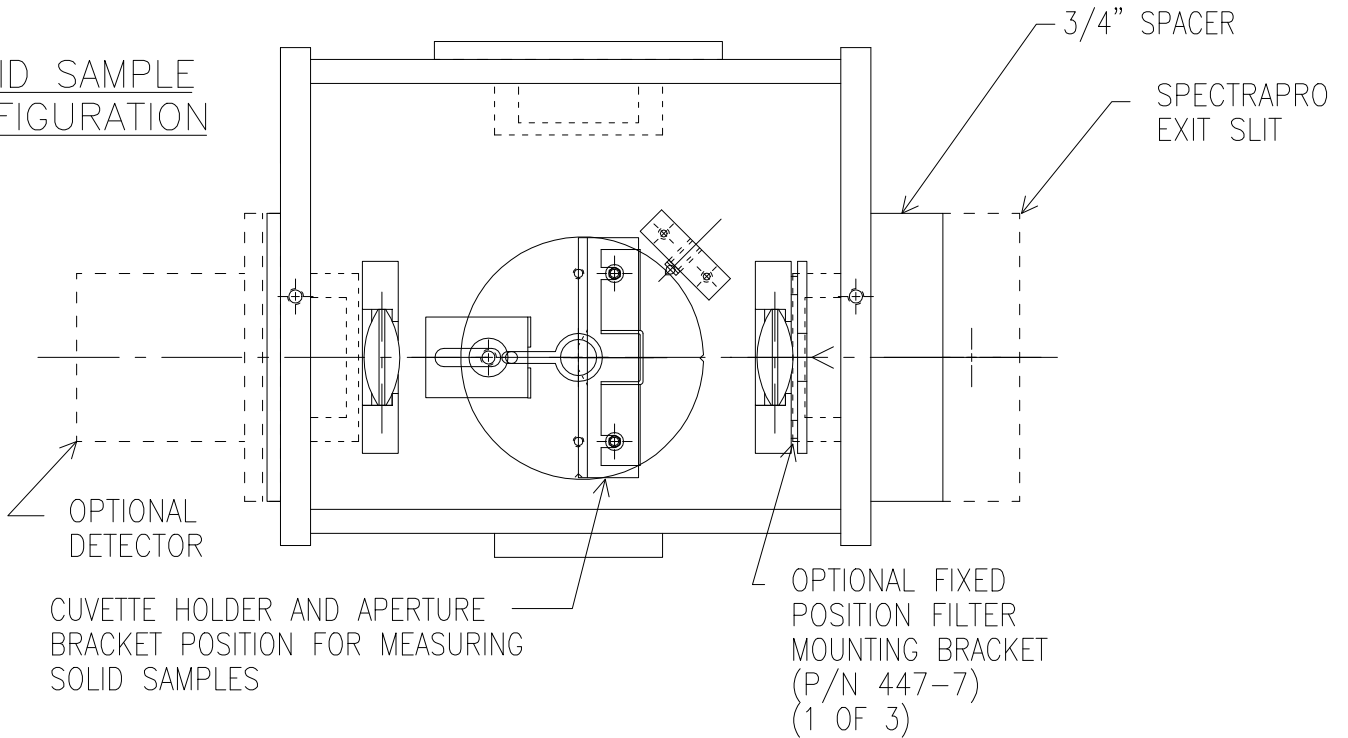

CUVETTE SAMPLE
 CONFIGURATION

